

Recombinant Anti-Cytokeratin 20 Rabbit mAb

Product Manual

Cat No: GMA1127

Source: Rabbit

Reactivity: H, M

Applications: WB, IH, IF/IC

***Application Key:**

H - Human, M - Mouse, R - Rat, B - Bovine, C - Chicken, D - Dog, G - Goat, Mk - Monkey, P - Pig, Rb - Rabbit, S - Sheep,
Z - Zebrafish

***Species Reactivity Key:**

E- ELISA, WB - Western blot, IH - Immunohistochemistry, IF - Immunofluorescence, FC - Flow cytometry, IC -
Immunocytochemistry, IP - Immunoprecipitation, ChIP - Chromatin Immunoprecipitation, EMSA - Electrophoretic
Mobility Shift Assay, BL - Blocking, SE - Sandwich ELISA, CBE - Cell-based ELISA, RNAi - RNA interference

Datasheet

Description: Recombinant rabbit monoclonal antibody to Cytokeratin 20

Immunogen: Recombinant fusion protein of human Cytokeratin 20. The exact sequence is proprietary.

Purification: The antibody was purified by immunogen affinity chromatography.

Clonality: Monoclonal

Form: Liquid in PBS, pH 7.3, 50% glycerol, and 0.05% Proclin300.

Dilution: WB (1/500 - 1/1000), IH (1/100 - 1/500), IF/IC (1/100 - 1/300)

Gene Symbol: KRT20


Alternative Names: Keratin, type I cytoskeletal 20; Cytokeratin-20; CK-20; Keratin-20; K20; Protein IT

Entrez Gene (Human): 54474;


SwissProt (Human): P35900;

Storage/Stability: Shipped at 4°C. Upon delivery aliquot and store at -20°C for one year. Avoid freeze/thaw cycles.


Validation

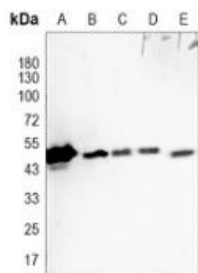

Western blot analysis of Cytokeratin 20 expression in HT29 (A), U87 (B), PC3 (C), mouse kidney (D), mouse spleen (E) whole cell lysates. (Predicted band size: 60 kD; Observed band size: 60 kD)

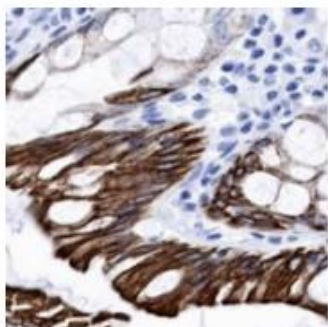

Immunohistochemical analysis of Cytokeratin 20 staining in human colon formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

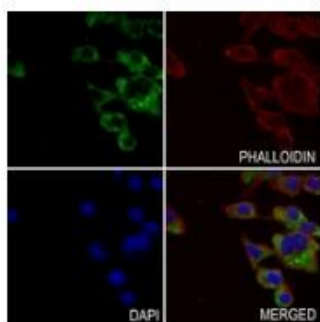

Immunofluorescent analysis of Cytokeratin 20 staining in HepG2 cells. Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at room temperature. Cells were probed with the primary antibody in 3% BSA-PBS and incubated overnight at 4 °C in a humidified chamber. Cells were washed with PBST and incubated with an AF488-conjugated secondary antibody (green) in PBS at room temperature in the dark. Phalloidin - AF594 was used to stain Actin filaments (red). DAPI was used to stain the cell nuclei (blue).

官方网址: <http://www.genescion.com.cn>

订货热线: 4006169114、020-84224925

Email: whiga22@126.com

